



RETRACTOR MAKER

PRODUCT BROCHURE OF

Bookwalter Hand Surgery Retractor System

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

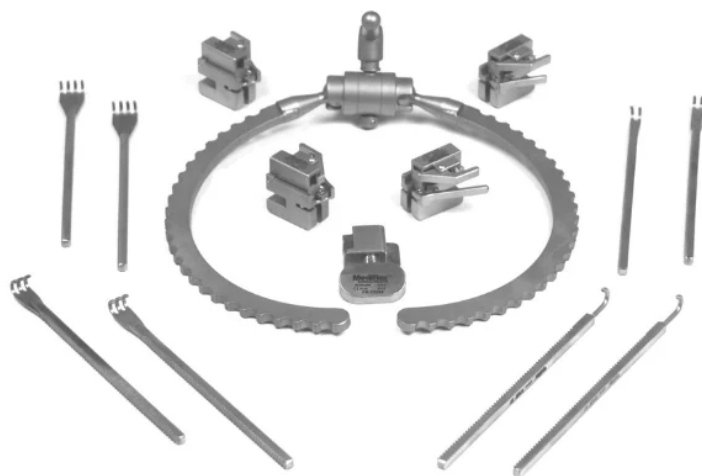

Bookwalter Hand Surgery Retractor System

SKU RM-BHSRS-1

\$999 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
2	Mini-Bookwalter Blade, Single Hook	Standard Mini — sharp or blunt tip
2	Mini-Bookwalter Blade, Double Hook	Standard Mini — sharp or blunt tip
2	Mini-Bookwalter Blade, Triple Hook	Standard Mini — sharp or blunt tip
2	Mini-Bookwalter Blade, Quad Hook	Standard Mini — sharp or blunt tip
1	Mini-Bookwalter Twist Ratchet	Standard Mini
4	Mini-Bookwalter Tilt Ratchet Mechanism	Standard Mini
1	Mini-Bookwalter Split Ring	7" (17.5 cm)
1	Sterilization Box	Protective Case

PRODUCT DESCRIPTION

The Bookwalter Hand Surgery Retractor System is a compact self-retaining retractor system designed to provide stable, hands-free retraction during delicate hand, wrist and upper extremity surgical procedures. The set is built around a 7" (17.5 cm) Mini-Bookwalter split ring, which forms the working frame around the operative field.

Four mini tilt ratchet mechanisms and a mini twist ratchet secure the included blades to the split ring and allow the angle and tension of each blade to be adjusted independently. The Bookwalter Hand Surgery Retractor System includes single hook, double hook, triple hook and quad hook Mini-Bookwalter blades, available with sharp or blunt tips, so the blade configuration can be selected to suit the exposure required. A StrongArm Holder is available separately for table mounting.

All components are manufactured from Premium Grade AISI 304/316 Stainless Steel (Non-Magnetic) and are supplied in a dedicated sterilization box for organized storage and processing. The system is reusable and fully autoclavable.

The Bookwalter Hand Surgery Retractor System is supplied non-sterile. It must be cleaned and sterilized according to the healthcare facility's validated procedures before use. Retractor Maker, based in Sialkot, Pakistan, manufactures and supplies this system directly to hospitals, surgeons and distributors.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractormaker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractormaker.com
retractormanufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)