



RETRACTOR MAKER

PRODUCT BROCHURE OF

Aufricht Nasal Retractor - Fiber Optic

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

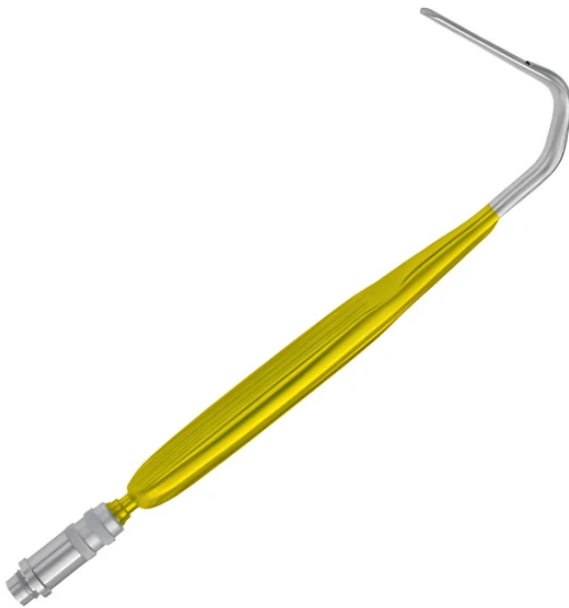

Aufricht Nasal Retractor - Fiber Optic

SKU RM-ANR-FO-1

\$149 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
1	Aufricht Nasal Retractor - Fiber Optic	6 ¼" (15.9 cm) length, 7.0 mm concave blade, light bundle
1	Universal Light Guide Adapter	Fits standard OR fiber-optic cables
1	Protective Sleeve / Box	Standard transit packaging

PRODUCT DESCRIPTION

The Aufricht Nasal Retractor - Fiber Optic is a premier illuminated rhinoplasty retractor engineered to provide crystal-clear, shadowless illumination and tissue elevation within deep intranasal corridors. By combining anatomical retraction with an integrated fiber-optic light channel, this instrument eliminates shadowing caused by overhead surgical lights, providing direct line-of-sight visualization to the osteotomy site and cartilaginous septum.

Designed to optimize working space in narrow nasal pathways, this Aufricht Nasal Retractor - Fiber Optic features an overall length of 6 ¼ inches (15.9 cm) and a slender 7.0 mm wide concave blade profile. The concave geometry creates additional clearance for osteotomes, rasps, and suction cannulas while lifting the dorsal flap. The instrument incorporates a seamless fiber-optic light guide terminated with a universal adapter compatible with standard hospital light cables (Storz, ACMI, Wolf).

Crafted from Premium Grade AISI 304/316 Stainless Steel (Non-Magnetic) with high-transmission optical fibers and an anti-glare satin finish, Retractor Maker produces the Aufricht Nasal Retractor - Fiber Optic non-sterile and fully autoclavable, delivering exceptional durability and brilliant illumination for demanding aesthetic nasal procedures.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractormaker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractormaker.com
retractormanufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)