



RETRACTOR MAKER

PRODUCT BROCHURE OF

Bookwalter Small Incision Retractor System

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

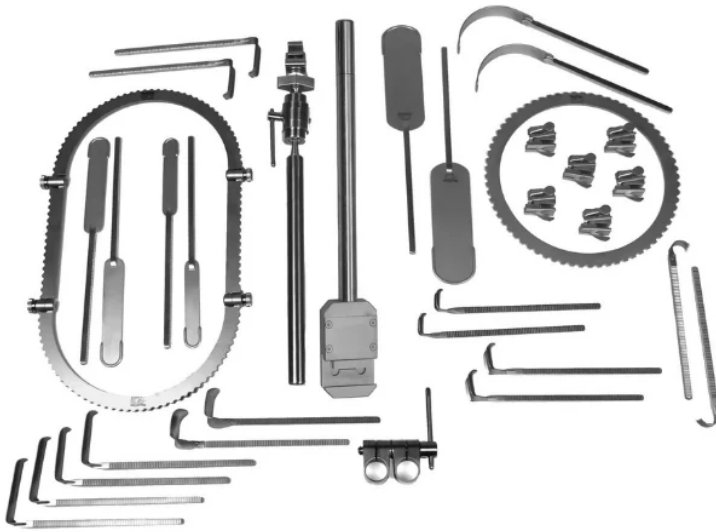

Bookwalter Small Incision Retractor System

SKU RM-BSIRS-1

\$2,199 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
1	Field Post	12" (30 cm) vertical support rod
1	Stainless Steel Post Coupler	1" × 1" (2.5 × 2.5 cm) heavy-duty rail clamp
1	Horizontal Flex Bar	1" Diameter (2.5 cm) multi-jointed arm
1	Small Round Solid Ring	8½" (21.3 cm) Internal Diameter (ID)
1	Small Oval Segmented Ring	8½" × 14" (20 × 35 cm) hinged reconfigurable frame
6	Tilt Ratchets	Low-profile mechanical frame locks
2	Malleable Blades	1" × 3½" (2.5 × 12.5 cm)
2	Malleable Blades	¾" × 4" (1.9 × 10 cm)
2	Malleable Blades	2" × 6" (5 × 15 cm)
2	Deaver Blades	1" × 3½" (2.5 × 8.8 cm) curved pelvic blades
2	Kelly Blades	2" × 1½" (5 × 3.8 cm) shallow right-angle blades
2	Kelly Blades	1" × 2½" (2.5 × 6.3 cm) right-angle blades
2	Kelly Blades	1" × 3" (2.5 × 7.5 cm) right-angle blades
2	Deep Rake Blades	¾" Deep × ⅔" Wide (1.9 × 1.7 cm) multi-prong claws

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
2	Richardson Blades	$\frac{3}{4}$ " $\times$ 2" (1.9 $\times$ 5 cm) shallow flat-lipped blades
2	Richardson Blades	1" $\times$ 1 $\frac{1}{4}$ " (2.5 $\times$ 3.1 cm) shallow flat-lipped blades
2	Richardson Blades	$\frac{3}{4}$ " $\times$ 1" (1.9 $\times$ 2.5 cm) shallow flat-lipped blades
1	Sterilization Box	Customized perforated organizer container

PRODUCT DESCRIPTION

The Bookwalter Small Incision Retractor System is a table-mounted, self-retaining retractor system. It anchors to the operating table and holds each blade in position, providing stable, hands-free surgical exposure during procedures performed through limited incisions where precise access is required while minimizing tissue disruption.

The mounting system comprises a 12" (30 cm) field post, a 1" $\times$ 1" (2.5 $\times$ 2.5 cm) heavy-duty stainless steel post coupler rail clamp and a 1" (2.5 cm) diameter multi-jointed horizontal flex bar. Two ring options are supplied: a small round solid ring with an 8 $\frac{1}{2}$ " (21.3 cm) internal diameter and a small oval segmented ring measuring 8 $\frac{1}{2}$ " $\times$ 14" (20 $\times$ 35 cm), which is hinged so the frame can be reconfigured around the incision.

Six low-profile tilt ratchets lock to the ring and hold each blade at the required depth and angle. The blade set includes malleable blades in three sizes, curved Deaver pelvic blades, shallow right-angle Kelly blades in three sizes, multi-prong deep rake blades and shallow flat-lipped Richardson blades in three sizes, allowing exposure to be built up progressively through a small opening.

All components are manufactured from Premium Grade AISI 304/316 Stainless Steel (Non-Magnetic). The Bookwalter Small Incision Retractor System is brand new, reusable and fully autoclavable, and is supplied non-sterile — cleaning and sterilization must follow the healthcare facility's validated procedures before use. The Bookwalter Small Incision Retractor System is supplied in a customized perforated sterilization box that organizes each component for processing and storage. Retractor Maker manufactures the system in Sialkot, Pakistan and supplies it directly to hospitals, surgeons and distributors, with a minimum order quantity of one full set.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractor maker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractor maker.com
retractor manufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)