



RETRACTOR MAKER

PRODUCT BROCHURE OF

Bookwalter Rectal/Trans-Anal Surgery Retractor

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

Donald Cameron

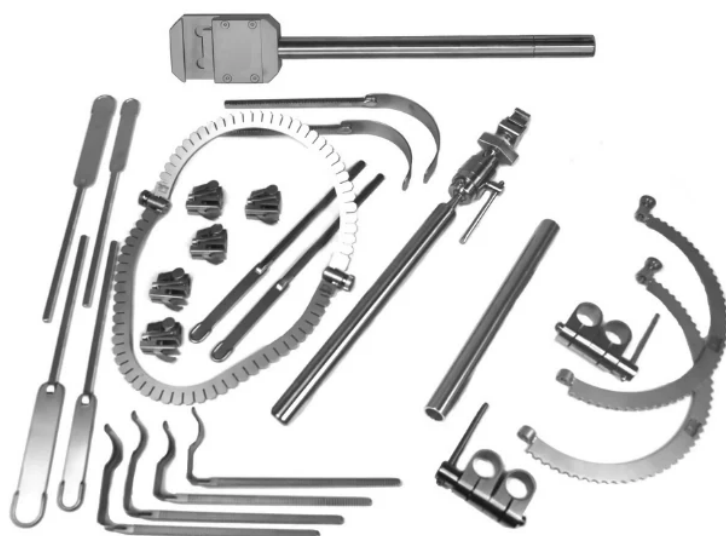

Bookwalter Rectal/Trans-Anal Surgery Retractor

SKU RM-BRTASR-1

\$2,199 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
1	Field Post	12" (30 cm)
2	Stainless Steel Post Couplers	1" x 1" (2.5 cm x 2.5 cm)
2	Horizontal Flex Bars	1" Diameter (2.5 cm)
1	Vertical Extension Bar	12" (30 cm)
5	Tilt Ratchets	Standard Multi-Angle
1	Small Half Circle Ring	8½" (21.3 cm) ID
1	Segmented Perineal Ring	9½" x 14½" (23.8 cm x 36.3 cm)
2	Rectal Blades	7⁄8" x 3⁷⁄8" (2.2 cm x 9.7 cm)
2	Rectal Blades	7⁄8" x 2⁷⁄8" (2.2 cm x 7.2 cm)
2	Malleable Blades	1" x 3½" (2.5 cm x 10 cm)
2	Malleable Blades	½" x 4" (1.3 cm x 8.8 cm)
2	Malleable Blades	¾" x 4" (1.9 cm x 10 cm)
2	Deaver Blades	1" x 3½" (2.5 cm x 8.8 cm)
1	Sterilization Box	Protective Case

PRODUCT DESCRIPTION

The Bookwalter Rectal/Trans-Anal Surgery Retractor is a table-mounted self-retaining retractor system designed to provide stable, hands-free surgical exposure during colorectal, trans-anal, rectal and perineal procedures. Its mounting framework uses a 12" (30 cm) field post, two 1" diameter horizontal flex bars, a 12" (30 cm) vertical extension bar and two 1" x 1" stainless steel post couplers, allowing flexible positioning of the frame relative to the operating table and patient position.

The working frame can be assembled from a 9½" x 14½" (23.8 x 36.3 cm) segmented perineal ring or an 8½" (21.3 cm) internal-diameter small half circle ring, depending on the approach. Five standard multi-angle tilt ratchets attach to the ring and allow each blade to be angled and tensioned independently.

The included blade set covers a range of exposure requirements: rectal blades in 7/8" x 37/8" and 7/8" x 27/8" sizes, malleable blades in 1" x 3½", ½" x 4" and ¾" x 4" sizes, and 1" x 3½" Deaver blades. All components are manufactured from Premium Grade AISI 304/316 Stainless Steel (Non-Magnetic) and are supplied in a dedicated sterilization box.

The Bookwalter Rectal/Trans-Anal Surgery Retractor is supplied non-sterile and is reusable and fully autoclavable. The Bookwalter Rectal/Trans-Anal Surgery Retractor must be cleaned and sterilized according to the healthcare facility's validated procedures before use. Retractor Maker manufactures the system in Sialkot, Pakistan and supplies it directly to hospitals, surgeons and distributors.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractormaker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractormaker.com
retractormanufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)