



RETRACTOR MAKER

PRODUCT BROCHURE OF

Bookwalter Posterior Spinal Retractor

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

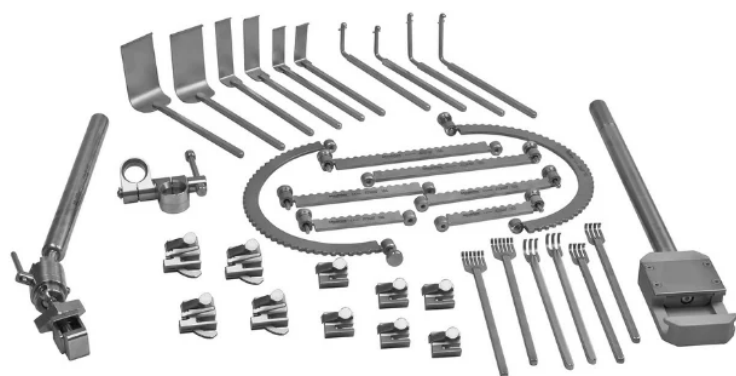

Bookwalter Posterior Spinal Retractor

SKU RM-BPSR-1

\$2,399 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
1	Field Post	12" (30.5 cm)
1	Post Coupler	Standard
1	Horizontal Bar	1" Diameter (2.5 cm)
6	Standard Ratchets	Non-Tilt
4	Tilt Ratchet Mechanisms	Multi-Angle
1	Half Circle Ring Segments	8½" Diameter (21.6 cm)
1	Straight Ring Segment	4" Diameter (10.2 cm)
1	Straight Ring Segment	6" Diameter (15.2 cm)
1	Straight Ring Segment	8" Diameter (20.3 cm)
2	Meyerding Blades	⅔" × 2⅓" (1.7 × 5.9 cm)
2	Meyerding Blades	1" × 3⅓" (2.5 × 8.5 cm)
2	Meyerding Blades	2" × 3¾" (5.1 × 9.5 cm)
2	Scoville Hooks	Small
2	Scoville Hooks	Large

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
2	3-Prong Rakes	$\frac{2}{3}$ " Wide $\times$ $\frac{3}{4}$ " Deep (1.7 $\times$ 1.9 cm)
2	4-Prong Rakes	$\frac{3}{4}$ " Wide $\times$ $\frac{1}{3}$ " Deep (1.9 $\times$ 0.8 cm)
2	6-Prong Rakes	1" Wide $\times$ $\frac{1}{3}$ " Deep (2.5 $\times$ 0.8 cm)
1	Sterilization Box	Protective Case

PRODUCT DESCRIPTION

The Bookwalter Posterior Spinal Retractor is a table-mounted self-retaining retractor set designed to provide stable, hands-free surgical exposure during posterior approaches to the cervical, thoracic and lumbar spine. The Bookwalter Posterior Spinal Retractor mounts to the operating table using a 12" (30.5 cm) field post, a standard post coupler and a 1" (2.5 cm) diameter horizontal bar.

The working frame is modular. An 8½" (21.6 cm) half circle ring segment combines with 4", 6" and 8" straight ring segments so the frame can be assembled to match the length and width of the incision. Six standard non-tilt ratchets and four multi-angle tilt ratchet mechanisms secure the blades to the frame and allow blade angle and retraction depth to be set independently.

The included retraction hardware covers a range of exposures: Meyerding blades in $\frac{2}{3}$ " $\times$ $2\frac{1}{3}$ ", 1" $\times$ $3\frac{1}{3}$ " and 2" $\times$ $3\frac{3}{4}$ " sizes, small and large Scoville hooks, and 3-prong, 4-prong and 6-prong rakes. All components are manufactured from Premium Grade AISI 304/316 Stainless Steel (Non-Magnetic) and supplied in a dedicated sterilization box.

The Bookwalter Posterior Spinal Retractor is supplied non-sterile and is reusable and fully autoclavable. It must be cleaned and sterilized according to the healthcare facility's validated procedures before use. Retractor Maker manufactures and supplies the set directly from Sialkot, Pakistan.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractormaker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractormaker.com
retractormanufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)