



RETRACTOR MAKER

PRODUCT BROCHURE OF

Bookwalter GYN Magrina Pelvic Surgery Retractor System

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

Donald Cameron

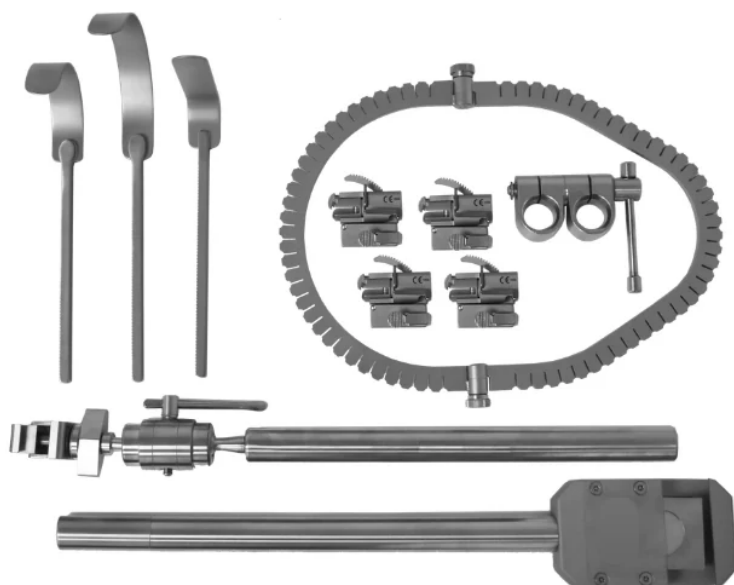

Bookwalter GYN Magrina Pelvic Surgery Retractor System

SKU RM-BGMPRS-1

\$1,499 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
1	Field Post	12" (30 cm)
1	Post Coupler	Standard
1	Horizontal Flex Bar	Standard Multi-Jointed
1	Perineal Ring with Notches	9½" × 14½" (24.1 × 36.8 cm)
4	Rotating Tilt Ratchet	Standard with Ring Retention
2	Vaginal Lateral Blades	1" × 3" (2.5 × 7.6 cm)
1	Deaver Blade	1" × 3½" (2.5 × 8.9 cm)
1	Deaver Blade	1" × 5" (2.5 × 12.7 cm)
1	Vaginal Posterior Blade	1" × 4½" (2.5 × 11.4 cm)
1	Sterilization Box	Protective Case

PRODUCT DESCRIPTION

The Bookwalter GYN Magrina Pelvic Surgery Retractor System is a table-mounted self-retaining retractor designed to provide stable, hands-free exposure during gynecological, vaginal and pelvic surgical procedures. Its modular framework combines a field post, horizontal flex bar, post coupler, notched perineal ring and rotating tilt ratchet mechanisms to support selected blades around the operative field.

The system is configured with a 9½" × 14½" notched perineal ring and four rotating tilt ratchets with ring retention. The included vaginal lateral, Deaver and vaginal posterior blades provide a range of blade configurations for surgical exposure requirements.

The complete system is supplied in a dedicated protective sterilization box for organized storage and processing. The Bookwalter GYN Magrina Pelvic Surgery Retractor System is supplied non-sterile and is designed to be reusable and autoclavable. Cleaning and sterilization should be performed according to the healthcare facility's validated procedures before use.

Retractor Maker supplies the Bookwalter GYN Magrina Pelvic Surgery Retractor System directly to hospitals, distributors and other B2B customers.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractormaker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractormaker.com
retractormanufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)