



RETRACTOR MAKER

PRODUCT BROCHURE OF

Logan Wound Spreader

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

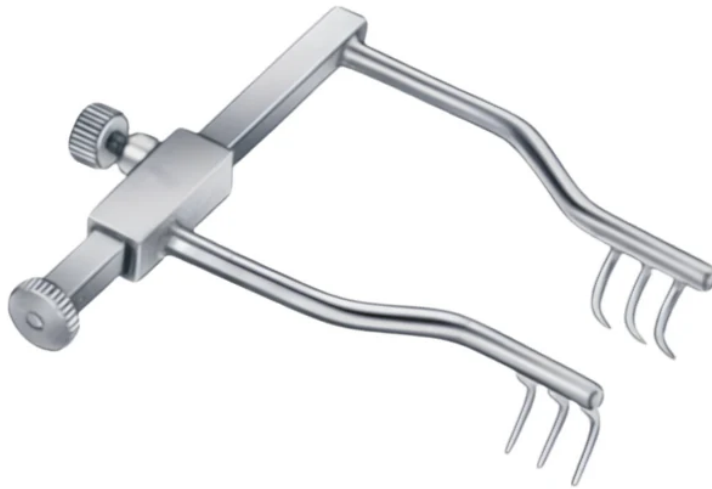

Logan Wound Spreader

SKU RM-LWS-4CM-1

\$199 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
1	Bottom Ratchet Guide Bar Frame	Main base track chassis serving as a guide for expanding tissue margins
1	Sliding Arm Bracket	Movable bracket component adjusting horizontally along the track
2	Screw Adjustment Mechanisms	Precision locking screw gears to shift and secure arms in position
2	Pronged Spreading Tips	Non-slip jaw ends featuring 3 × 3 sharp outward-curved teeth
1	Protective Storage / Box	Standard transit packaging

PRODUCT DESCRIPTION

The Logan Wound Spreader is a specialized compact self-retaining surgical retractor engineered specifically to provide secure, hands-free soft-tissue and wound margin exposure during ophthalmic, facial plastic, and micro-reconstructive surgery. Designed with an ultra-short horizontal ratchet bar, this instrument allows surgeons to spread shallow incisions with calibrated micro-precision, eliminating assistant hand placement while maintaining an unobstructed field under operating loupes and microscopes.

Built specifically for delicate anatomical corridors, this Logan Wound Spreader configuration features an overall instrument length of 4 cm (1 ½ inches, \$199). Handcrafted from Premium Grade AISI 304/316 Stainless Steel (Non-Magnetic), the working end configuration incorporates non-slip jaw ends armed with 3 × 3 sharp outward-curved teeth specifically designed to safely anchor into dense aponeurosis, periosteum, and muscle margins without slipping. The sliding arm bracket shifts smoothly along the bottom ratchet guide bar frame and

locks rigidly via two precision screw adjustment mechanisms to eliminate user fatigue.

Forged with an anti-glare satin matte finish, Retractor Maker manufactures the Logan Wound Spreader to withstand repetitive hospital autoclave cycles, providing ophthalmic surgery centers, plastic clinics, and cranial microsurgery suites worldwide with exceptional structural rigidity, corrosion resistance, and complete autoclave durability.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractormaker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractormaker.com
retractormanufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)