



RETRACTOR MAKER

PRODUCT BROCHURE OF

Tupper's Universal Hand Retractor

CERTIFIED & APPROVED

✓ EN ISO 13485:2016 + A11:2021



ISO 9001:2015



EN ISO 7153-1:2016



cGMP Approved




SCCI Registered since 2023



SIMAP Approved since 2023

All procedures and processes follow FDA guidelines.

Donald Cameron

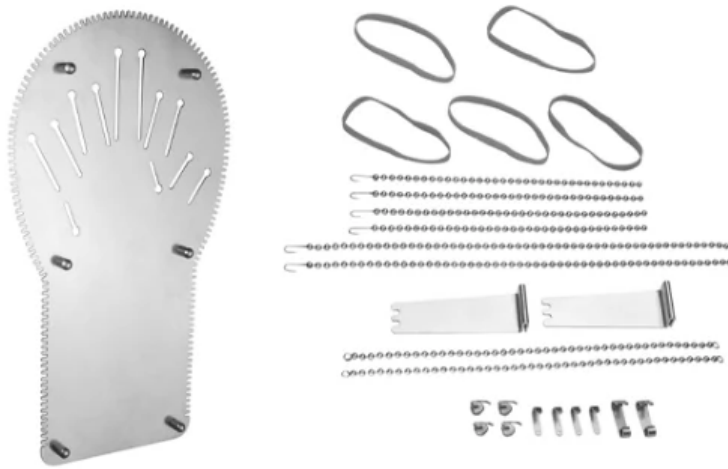

Tupper's Universal Hand Retractor

SKU RM-TUHR-1

\$699 USD

COMPONENT LIST

QUANTITY	COMPONENT	DIMENSIONS / DESCRIPTION
1	Hand Table	Central operative surface
2	Elevated Chain Holders	Frame perimeter supports
2	Chain Holder Hooks	Fixed attachment points
2	Sliding Chain Holder Hooks	Adjustable attachment points
4	Ball Chains with Skin Hooks	21 cm
6	Ball Chains for Arm/Wrist Fixation	26 cm
1 Set	Rubber Bands	5 pieces

PRODUCT DESCRIPTION

Tupper's Universal Hand Retractor is an internationally acclaimed surgical retractor platform engineered to provide precise, hands-free positioning and multi-directional tissue retraction during delicate hand and microsurgical procedures. By anchoring the patient's hand and digits to a central operative table via elastic and ball-chain tension, this system frees surgical assistants while establishing an unencumbered microsurgical field.

Designed for complete procedural adaptability, this Tupper's Universal Hand Retractor set features a solid stainless steel hand table, two elevated perimeter chain holders, two fixed chain holder hooks, and two sliding chain holder hooks that lock anywhere along the frame. The set includes four 21 cm ball chains fitted with sharp skin hooks, six 26 cm ball chains for wrist and arm immobilization, and a 5-piece set of medical-grade rubber bands for elastic digital counter-traction.

Manufactured entirely from Premium Grade AISI 304/316 Stainless Steel (Non-Magnetic) with a reflection-free satin finish, Retractor Maker supplies Tupper's Universal Hand Retractor non-sterile and fully autoclavable inside an individual protective packaging box, delivering exceptional clinical durability for hand trauma and reconstructive centers worldwide.



THANK YOU

We look forward to discussing your surgical retractor requirements.

CONTACT RETRACTOR MAKER

EMAIL

info@retractormaker.com

WHATSAPP

+92 327 8345101

WEBSITES

retractormaker.com
retractormanufacturer.com

ADDRESS

Retractor Maker, Near FFRV+F29, Roras Road, Sialkot-51310, Pakistan

FOLLOW US



[LinkedIn](#)



[Facebook](#)



[Instagram](#)



[X \(Twitter\)](#)



[YouTube](#)



[Pinterest](#)